



Genomic DNA from plant

User manual

NucleoSpin[®] 96 Plant II

NucleoSpin[®] 96 Plant II Core Kit

August 2020 / Rev. 04

Table of contents

1 Components	4
1.1 Kit contents	4
1.2 Reagents to be supplied by user	5
2 Product description	6
2.1 The basic principle	6
2.2 Kit specifications	6
2.3 Required hardware	7
2.4 Recommended accessories for use of the NucleoSpin® 96 Plant II Core Kit	8
2.5 Automated processing on robotic platforms	10
2.6 Storage and homogenization of samples	10
2.7 Elution procedures	11
2 Product description	6
3 Storage conditions and preparation of working solutions	13
4 Safety instructions	14
5 Protocols	16
5.1 NucleoSpin® 96 Plant II – centrifuge processing	16
5.2 NucleoSpin® 96 Plant II – vacuum processing	20
6 Appendix	26
6.1 Troubleshooting	26
6.2 Ordering information	28
6.3 Product use restriction / warranty	29

1 Components

1.1 Kit contents

NucleoSpin® 96 Plant II			
	2 x 96 preps	4 x 96 preps	24 x 96 preps
REF	740663.2	740663.4	740663.24 ¹
Lysis Buffer PL1	125 mL	250 mL	6 x 250 mL
Lysis Buffer PL2 ²	100 mL	200 mL	6 x 200 mL
Precipitation Buffer PL3	25 mL	50 mL	6 x 50 mL
Binding Buffer PC	125 mL	250 mL	6 x 250 mL
Wash Buffer PW1	100 mL	2 x 100 mL	12 x 100 mL
Wash Buffer PW2 (Concentrate) ²	100 mL	2 x 100 mL	12 x 100 mL
Elution Buffer PE ³	60 mL	125 mL	6 x 125 mL
RNase A (lyophilized) ²	30 mg	2 x 30 mg	12 x 30 mg
NucleoSpin® Plant II Binding Plate (dark green rings)	2	4	24
MN Wash Plate	2	4	24
Rack of Tube Strips ⁴ (for lysis and elution)	4	8	32
Cap Strips	24	48	288
MN Square-well Block	6	12	72
Gas-permeable Foil	10	20	120
User manual	1	1	1

¹ The kit for 24 x 96 preparations REF 740663.24 consists of 6 x REF 740663.4.

² For preparation of working solutions and storage conditions see section 3.

³ Composition of Elution Buffer PE: 5 mM Tris/HCl, pH 8.5

⁴ 1 rack = 12 strips with 8 tubes each, Cap Strips included

Kit contents continued

NucleoSpin® 96 Plant II Core Kit	
REF	4 x 96 preps 740468.4
Lysis Buffer PL1	250 mL
Lysis Buffer PL2 ¹	200 mL
Precipitation Buffer PL3	50 mL
Binding Buffer PC	250 mL
Wash Buffer PW1	2 x 100 mL
Wash Buffer PW2 (Concentrate) ¹	2 x 100 mL
Elution Buffer PE ²	125 mL
RNase A (lyophilized) ¹	2 x 30 mg
NucleoSpin® Plant II Binding Plate (dark green rings)	4
User manual	1

1.2 Reagents to be supplied by user

- 96–100 % ethanol

¹ For preparation of working solutions and storage conditions see section 3.

² Composition of Elution Buffer PE: 5 mM Tris/HCl, pH 8.5

2 Product description

2.1 The basic principle

The **NucleoSpin® 96 Plant II** kit is designed for the isolation of genomic DNA from plant materials. After the plant samples have been homogenized, the DNA can be extracted with lysis buffers containing chaotropic salts, denaturing agents and detergents. The standard isolation ensures the lysis of plant material with the CTAB Lysis Buffer PL1, which is specially developed for plants. In addition, an SDS based lysis buffer, Buffer PL2, is provided as an alternative. Buffer PL2 requires subsequent protein precipitation with potassium acetate. Lysates should be cleared by centrifugation in order to remove polysaccharides, contaminations, and residual cellular debris. The clear supernatant is mixed with Binding Buffer PC to create conditions for optimal binding to the silica membrane in the binding plate. After washing with two different buffers (Buffer PW1 and Buffer PW2) DNA can be eluted in low salt Buffer PE or water and is ready to use for subsequent analysis and processing.

2.2 Kit specifications

- **NucleoSpin® 96 Plant II** is designed for the isolation of genomic DNA from plant material.
- **NucleoSpin® 96 Plant II** allows parallel purification of multiples of 96 samples each with up to 100 mg sample per well (wet weight).
- Depending on the individual sample, **NucleoSpin® 96 Plant II** shows yields in the range of 1–30 µg DNA (maximum column capacity is about 30 µg) with an A_{260}/A_{280} ratio between 1.80 and 1.90 and typical concentrations of 100–200 ng/µL. The amount of DNA that can be expected per mg of sample extracted depends on the size and ploidy of the genome. For example, 100 mg fresh wheat with a hexaploid genome (1.7 x 10¹⁰ bp) contain 30 µg DNA, whereas the same amount of Arabidopsis with a smaller diploid genome (1.9 x 10⁸ bp) yields only 3 µg DNA.
- The eluted DNA is ready to use in subsequent reactions like PCR, restriction analysis, etc.
- **NucleoSpin® 96 Plant II** can be processed under vacuum or in a centrifuge.
- Two lysis buffers, based on CTAB (PL1) or SDS (PL2) are provided.
- **NucleoSpin® 96 Plant II** can be used manually with the NucleoVac 96 Vacuum Manifold (see ordering information) or other vacuum devices.

Kit specifications at a glance

Parameter	NucleoSpin® 96 Plant II
Technology	Silica membrane technology
Format	96-well plates
Processing	Manual or automated, vacuum or centrifugation
Sample material	20–100 mg plant tissue, plant cells (wet weight)

Kit specifications at a glance

Fragment size	50 bp–approx. 50 kpb
Typical yield	1–30 µg
A ₂₆₀ /A ₂₈₀	1.8–1.9
Elution volume	100–200 µL
Preparation time	60 min/plate (excl. lysis)
Binding capacity	30 µg

2.3 Required hardware

NucleoSpin® 96 Plant II can be processed under vacuum or with centrifugation. Certain hardware for processing is required.

Centrifugation

For centrifugation, a microtiterplate centrifuge is required. This centrifuge must be able to accommodate the NucleoSpin® Plant II Binding Plate stacked on a Round- or Square-well Block and reach accelerations of 5,600–6,000 x g is required (bucket height: 85 mm).

Regarding waste collection, suitable consumables (e.g., MN Square-well Blocks) are necessary and they are not included in the kit. For the most convenient handling, without the need of emptying and reusing the MN Square-well Blocks, we recommend using six MN Square-well Blocks if two 96-well plates are processed at once (see ordering information). Alternatively, it is possible to empty the MN Square-well Blocks after every centrifugation step, reducing the amount of MN Square-well Blocks needed.

Vacuum processing

The **NucleoSpin® 96 Plant II** kit can be used with the NucleoVac 96 Vacuum Manifold (see ordering information). When using **NucleoSpin® 96 Plant II** with less than 96 samples, Self adhering PE Foil (see ordering information) should be used in order to close and protect non-used wells of the NucleoSpin® Plant II Binding Plate and thus guarantee proper vacuum.

Establish a reliable vacuum source for the NucleoVac 96 Vacuum Manifold. The manifold may be used with a vacuum pump, house vacuum, or water aspirator. We recommend a vacuum of -0.2 to -0.4 bar (reduction of atmospheric pressure). The use of the NucleoVac Vacuum Regulator (see ordering information) is recommended. Alternatively, adjust the vacuum so that during the purification the sample flows through the column with a rate of 1–2 drops per second. Depending on the amount of sample being used, the vacuum times may need to be increased for complete filtration.

Additionally, a suitable centrifuge for sample preparation steps may be required.

For general consumables and equipment needed, please see section 1.2.

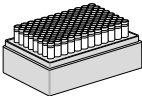
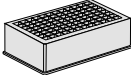
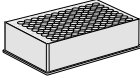
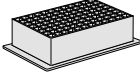
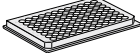
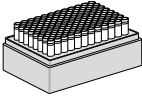
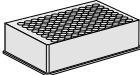
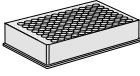
2.4 Recommended accessories for use of the NucleoSpin® 96 Plant II Core Kit

The **NucleoSpin® 96 Plant II Core Kit** provides all necessary buffers, enzymes and NucleoSpin® Binding Plates. Accessories (e.g., lysis plates, waste collection plates, elution plates or tubes) are not provided with the Core Kit. The reduced kit composition along with a large variety of separately available accessories, allow optimal adjustment of the kit to individual user needs. The user can select additional consumables according to his requirements for highest flexibility.

The **NucleoSpin® 96 Plant II Core Kit** provides buffers RNase A, and NucleoSpin® Binding Plates only. Accessory plates (e.g., elution plates) are not provided with the core kit. The user can individually select additional consumables from a variety of suitable accessory plates according to his requirements for highest flexibility.

For use of **NucleoSpin® 96 Plant II Core Kit** follow the standard protocols (see section 5.1 or 5.2, respectively).

Recommended accessories for use of the **NucleoSpin® 96 Plant II Core Kit** are available from MACHEREY-NAGEL. For ordering information please refer to section 6.2.

	Suitable consumables, not supplied with the core kits		Remarks
1. Homogenize samples	Rack of Tube Strips with Cap Strips		
4. Adjust binding conditions	Square-well Block or		For mixing cleared lysate with Buffer PC
	Round-well Block or		
	MN Square-well Block		
7. Wash silica membrane	MN Wash Plates		MN Wash Plate minimizes the risk of cross contamination (vacuum processing)
8. Elute DNA	Rack of Tubes Strips with Cap Strips or		
	Round-well Block or		
	Round-well Block Low (centrifugation only)		

2.5 Automated processing on robotic platforms

NucleoSpin® 96 Plant II can be used fully automated on many common laboratory workstations. For the availability of scripts and general considerations about adapting **NucleoSpin® 96 Plant II** on a certain workstation please contact MN. Full processing under vacuum enables complete automation without the need for centrifugation steps for drying of the binding membrane or for elution.

The risk of cross-contamination is reduced by optimized vacuum settings during the elution step and by the improved shape of the outlets of the NucleoSpin® Plant II Binding Plate.

Drying of the NucleoSpin® Plant II Binding Plates under vacuum is sufficient because the bottom of the plate is protected by the MN Wash Plate during the washing steps. As a result, it is recommended to integrate the MN Wash Plate into the automated procedure to protect against these wash buffer residues. The MN Frame (see ordering information) can be used to position the disposable MN Wash Plate inside the vacuum chamber. This also reduces the risk of cross-contamination as common metal adaptors tend to get contaminated by gDNA. In addition, thorough cleaning of the vacuum chamber is recommended after each run to prevent forming of gDNA-containing aerosols.

Visit MN online at www.mn-net.com or contact your local MACHEREY-NAGEL distributor for technical support regarding hardware, software, setup instructions, and selection of the protocol. Several application notes of the **NucleoSpin® 96 Plant II** kit on various liquid handling instruments can also be found at www.mn-net.com at Bioanalysis / Literature.

2.6 Storage and homogenization of samples

We recommend using young plant samples and keeping the plants in the dark for about 12 h before collecting samples (if possible) in order to reduce the polysaccharide content.

Plant samples can be stored frozen, under ethanol, or lyophilized. In many cases lyophilized, dried material can be processed more easily and gives higher yield. However, keep in mind that dried samples may reduce the amount of starting material by the factor 5 (for example, 20 mg dried plant leaves vs. 100 mg fresh weight).

As plant tissue is very robust, the lysis procedure is most effective with well homogenized, powdered samples. Suitable methods include grinding with pestle and mortar in the presence of liquid nitrogen or using steel beads. We also recommend the use of other commercial homogenizers, bead mills, etc.

Methods to homogenize samples

- Commercial homogenizers, for example Crush Express for 96-well homogenization (Saaten-Union Resistenzlabor GmbH, D-33818 Leopoldshöhe), Tissue Striker (www.KisanBiotech.com), or Geno/Grinder 2000 can be used.
- Samples can be disrupted using bead based homogenization tools, for example, GenoGrinder (<http://www.spexcsp.com> or for Germany www.c3-analysentechnik.de) or Mixer Mill MM400 (<http://www.retsch.com/products/milling/ball-mills/mm-400/>). Please refer to instrument manufacturers recommendations for suitable plates or tubes for homogenization.
- Homogenizing samples by VA steel beads (diameter: 3 mm): Put 4–5 beads and plant material together into a 15 ml plastic tube (Falcon), chill the tube in liquid nitrogen, and vortex for about 30 seconds (e.g., with a Multi Pulse Vortexer, Schütt Labortechnik GmbH, Postfach 3454, D-37024 Göttingen, Germany). Repeat this chilling and

vortexing procedure until the entire plant material is ground to a powder. Chill the tube once more and remove the beads by rolling them out gently or remove them with a magnet. Keep the material frozen throughout the whole homogenization procedure. Do not add nitrogen to the tube! This leads to sticking and loss of plant material attached to the beads. Alternatively, prefilled and ready to use MN Bead Tubes Type G (5 mm steel beads) allow an easy and convenient sample homogenization procedure.

- High throughput homogenization: Add the plant tissue to the individual tubes of the Tube Strips. Add one 3 mm stainless steel bead to each tube and close the individual tubes with Cap Strips. Freeze the sample in liquid nitrogen and insert the Rack of Tube Strips in a suitable homogenization tool (e.g., mixer mill). For disruption, shake the samples for 60–90 s at 30 Hz or until a homogenous plant powder has been formed. If necessary, repeat shaking once. Fresh plant material can also be homogenized with lysis buffer, however, homogenization of fresh plant material with lysis buffer may cause shearing of DNA. For frozen plant material thawing should be avoided during the homogenization. Samples should be frozen in liquid nitrogen before homogenization. Lyophilized or silica-gel dried material can be homogenized with or without lysis buffer. Homogenization of lyophilized tissue with lysis buffer may result in higher yield but also may cause shearing of DNA. Alternatively, prefilled and ready to use MN Bead Plate Type D allows an easy and convenient sample homogenization procedure.

2.7 Elution procedures

It is possible to adjust the elution method and the volume of the elution buffer to the specific application of interest. In addition to the standard method (recovery rate about 80–90 %) described in the protocols, there are 3 modifications possible:

- High yields: 90–100 % of bound nucleic acids can be eluted by performing two elution steps with volumes as indicated in the protocol, for example 2 x 100 µL. Finally, combine eluates and measure yield.
- Alternatively, use preheated Elution Buffer PE (70 °C): Preheat elution buffer to increase yield. After loading half of the preheated elution buffer (50 µL) onto the membrane, incubate the NucleoSpin® Plant II Binding Plate for 3 min at 60–70 °C. Centrifuge for elution as indicated. Repeat the elution step once.
- Highly concentrated eluates: Using a minimal elution volume (about 50 µL) about 70–80 % of bound nucleic acids can be eluted, resulting in highly concentrated eluates.

Elution may also be performed with Tris-EDTA-buffer (TE) with a pH equal or higher than 8. This will increase DNA stability during long term or multi-use storage at 4 °C (or ambient temperature) by inhibiting omnipresent DNases. However, EDTA interferes, depending on the final concentration, with certain downstream applications.

For optimal performance of isolated DNA in downstream applications, we recommend eluting with the supplied elution buffer and storing it, especially long term, at -20 °C. Several freeze-thaw cycles will not interfere with most downstream applications.

Performance of long-range PCR (e.g., > 10 kbp), or the detection limit of trace amounts of DNA species, may be reduced after multiple freeze-thaw cycles or prolonged storage of eluted DNA at 4 °C or room temperature. This is due to shearing of DNA or adsorption to surfaces.

Due to the dead volume of the silica membrane please note that the difference between the dispensed elution volume and the recovered elution buffer is approximately 45 µL (recovered elution volume = dispensed elution volume - 45 µL).

3 Storage conditions and preparation of working solutions

Attention: Buffer PL1 contains CTAB, Buffer PL2 contains SDS, Buffers PC and PW1 contain chaotropic salt! Wear gloves and goggles!

CAUTION: Buffers PC and PW1 contain guanidine hydrochloride which can form highly reactive compounds when combined with bleach (sodium hypochlorite). DO NOT add bleach or acidic solutions directly to the sample-preparation waste.

- Store **RNase A** at 4 °C on arrival (storage at 4 °C may cause precipitation of salts in different buffers).
- All other components can be stored at room temperature (18–25 °C) and are stable for up to 1 year.

Before starting any **NucleoSpin® 96 Plant II** protocol prepare the following:

- **Lysis Buffer PL2:** Check for precipitated SDS especially after storage at temperatures below 20 °C. If necessary incubate the bottle for several minutes at 30–40 °C and mix well until the precipitate is redissolved completely.
- **Wash Buffer PW2:** Add the indicated volume of ethanol (96–100 %) to Buffer PW2 Concentrate before first use. Store Buffer PW2 at room temperature (18–25 °C) for up to one year.
- **RNase A:** Add the given volume of water (indicated on the vial, see below) to lyophilized RNase A. Store the RNase A solution at 4 °C for up to 3 months. For longer storage (up to 1 year), the RNase A solution should be divided into small aliquots and stored at -20 °C.

	NucleoSpin® 96 Plant II	NucleoSpin® 96 Plant II	NucleoSpin® 96 Plant II	NucleoSpin® 96 Plant II Core Kit
	2 x 96 preps	4 x 96 preps	24 x 96 preps	4 x 96 preps
REF	740663.2	740663.4	740663.24	740468.4
Wash Buffer PW2 (Concentrate)	100 mL Add 400 mL ethanol	2 x 100 mL Add 400 mL ethanol to each bottle	12 x 100 mL Add 400 mL ethanol to each bottle	2 x 100 mL Add 400 mL ethanol to each bottle
RNase A (lyophilized)	30 mg Add 2.5 mL H ₂ O	2 x 30 mg Add 2.5 mL H ₂ O to each vial	12 x 30 mg Add 2.5 mL H ₂ O to each vial	2 x 30 mg Add 2.5 mL H ₂ O to each vial

4 Safety instructions

The following components of the **NucleoSpin® 96 Plant II** and **NucleoSpin® 96 Plant II Core** kits contain hazardous contents.

Wear gloves and goggles and follow the safety instructions given in this section.

Only harmful features need not be labeled with H and P phrases until 125 mL or 125 g.
Mindergefährliche Eigenschaften müssen bis 125 mL oder 125 g nicht mit H- und P-Sätzen gekennzeichnet werden.

Component	Hazard contents	GHS symbol	Hazard phrases	Precaution phrases
Inhalt	Gefahrstoff	GHS-Symbol	H-Sätze	P-Sätze
PC	guanidine hydrochloride 24–36 % and ethanol 35–55 % <i>Guanidinthiocyanat und Ethanol 35–55 %</i> CAS 50-01-, 64-17-5	 WARNING ACHTUNG	226, 302	210, 264W, 301+312, 330
PW1	guanidine hydrochloride 36–50 % and 2-propanol 20–35 % <i>Guanidinhydrochlorid 36–55 % und 2-Propanol 20–35 %</i> CAS 50-01-1, 67-63-0	 WARNING ACHTUNG	226, 302, 319, 336	210, 260D, 264W, 280sh, 301+312, 330
RNase A	RNase 90–100 % <i>RNase 90–1500 %</i> CAS 9001-99-4	 DANGER GEFAHR	334	261sh, 342+311

Hazard phrases

H 226	Flammable liquid and vapor. <i>Flüssigkeit und Dampf entzündbar.</i>
H 302	Harmful if swallowed. <i>Gesundheitsschädlich bei Verschlucken.</i>
H 319	Causes serious eye irritation. <i>Verursacht schwere Augenreizung.</i>
H 334	May cause allergy or asthma symptoms or breathing difficulties if inhaled. <i>Kann bei Einatmen Allergie, asthmaartige Symptome oder Atembeschwerden verursachen.</i>
H 336	May cause drowsiness or dizziness. <i>Kann Schläfrigkeit und Benommenheit verursachen.</i>

Precaution phrases

P 210	Keep away from heat/sparks/open flames/hot surfaces. No smoking. <i>Von Hitze, heißen Oberflächen, Funken, offenen Flammen sowie anderen Zündquellenarten fernhalten. Nicht rauchen.</i>
-------	---

P 260D	Do not breathe vapors. <i>Dampf nicht einatmen.</i>
P 261sh	Avoid breathing dust/vapors. <i>Einatmen von Staub/Dampf vermeiden</i>
P 264W	Wash with water thoroughly after handling. <i>Nach Gebrauch mit Wasser gründlich waschen.</i>
P 280sh	Wear protective gloves/eye protection. <i>Schutzhandschuhe/Augenschutz tragen.</i>
P 301+312	IF SWALLOWED: Call a POISON CENTER/doctor if you feel unwell. <i>BEI VERSCHLUCKEN: Bei Unwohlsein GIFTINFORMATIONSZENTRUM/Arzt anrufen.</i>
P 330	Rinse mouth. <i>Mund ausspülen.</i>
P 342+311	If experiencing respiratory symptoms: Call a POISON CENTER/doctor. <i>Bei Symptomen der Atemwege: GIFTINFORMATIONSZENTRUM/Arzt anrufen.</i>



The symbol shown on labels refers to further safety information in this section.
Das auf Etiketten dargestellte Symbol weist auf weitere Sicherheitsinformationen dieses Kapitels hin.

For further information please see Material Safety Data Sheets (www.mn-net.com).
Weiterführende Informationen finden Sie in den Sicherheitsdatenblättern (www.mn-net.com).

5 Protocols

5.1 NucleoSpin® 96 Plant II – centrifuge processing

- For hardware requirements, refer to section 2.3
- For detailed information on each step, see page 17.
- For use of the NucleoSpin® 96 Plant II Core Kit (REF 740468.4), refer to section 2.4 regarding recommended accessories.

Before starting the preparation:

- Check if Buffer PW2 and RNase A were prepared according to section 3.
- Set incubator or oven to 65 °C.
- Equilibrate Buffer PE to 70 °C.

Protocol at a glance

1	Homogenize samples	Up to 100 mg wet or 20 mg lyophilized plant tissue 5,600–6,000 x g, 2 min
2a	Cell lysis using Buffer PL1	500 µL PL1 10 µL RNase A Mix 65 °C, 30 min Proceed with step 3
2b	Cell lysis using Buffer PL2 and PL3	400 µL PL2 10 µL RNase A Mix 65 °C, 30 min 100 µL PL3 Mix and incubate on ice for 5 min Proceed with step 3
3	Clear lysate by centrifugation	5,600–6,000 x g, 20 min
4	Adjust binding conditions	Mix 450 µL PC with 400 µL cleared lysate
5	Transfer lysate to NucleoSpin® Plant II Binding Plate	

6	Bind DNA to silica membrane of the NucleoSpin® Plant II Binding Plate	5,600–6,000 x g, 2 min
7	Wash and dry silica membrane	400 µL PW1 5,600–6,000 x g, 2 min 700 µL PW2 5,600–6,000 x g, 2 min 700 µL PW2 5,600–6,000 x g, 10 min
8	Elute DNA	100 µL PE (70 °C) (incubate 2 min) 5,600–6,000 x g, 2 min Repeat once

Detailed protocol

- For hardware requirements, refer to section 2.3.
- For use of the NucleoSpin® 96 Plant II Core Kit (REF 740468.4), refer to section 2.4 regarding recommended accessories.

Before starting the preparation:

- Check if Buffer PW2 and RNase A were prepared according to section 3.
- Set incubator or oven to 65 °C.
- Equilibrate Buffer PE to 70 °C.

1 Homogenize samples

Fill **up to 100 mg wet plant tissue** (or up to **20 mg dried material**, for example lyophilized plant tissue) into each tube of the Tube Strips. Add one 3 mm diameter steel bead to each tube. Close the tubes with Cap Strips. Freeze samples in liquid nitrogen. Disrupt cells by vigorous shaking using a mixer mill. Centrifuge at **5,600 x g** for **2 min** and remove Cap Strips.

For further processing use either Buffer PL1 (2 a) or Buffers PL2 / PL3 (2 b)!

2a Cell lysis using Buffer PL1

Add **500 µL Buffer PL1** and **10 µL RNase A** to each sample. Close tubes again using new Cap Strips (supplied). Mix by vigorous **shaking** for **15–30 s**. **Spin briefly** for **30 s** at **1,500 x g** to collect any sample from the Cap Strips. Incubate samples at **65 °C** for **30 min**.

Depending on plant sample and available methods, Buffer PL1 and RNase A may be added to the plant material before homogenization by the appropriate mechanical method.

Proceed with step 3.

2b Cell lysis using Buffer PL2 and PL3

Add **400 µL Buffer PL2** and **10 µL RNase A** to each sample. Close tubes again using new Cap Strips (supplied). Mix by vigorous **shaking** for **15–30 s**. Spin briefly for **30 s** at **1,500 x g** to collect any sample from the Cap Strips. Incubate samples at **65 °C** for **30 min**.

Depending on plant sample and available methods, Buffer PL2 and RNase A may be added to the plant material before homogenization by the appropriate mechanical method.

Open tubes, add **100 µL Buffer PL3**, close tubes, mix thoroughly, and incubate for 5 min on ice to precipitate SDS completely.

3 Clear lysate by centrifugation

Centrifuge the samples for **20 min** at **full speed (5,600–6,000 x g)**. Remove Cap Strips.

4 Adjust binding conditions

Predispense **450 µL Binding Buffer PC** to each well of a MN Square-well Block. Add **400 µL cleared lysate** of each sample and mix by repeated pipetting up and down. Mix at least 3 times.

5 Transfer lysate to NucleoSpin® Plant II Binding Plate

Place NucleoSpin® Plant II Binding Plate on a MN Square-well Block. Transfer samples from the previous step into the wells of the NucleoSpin® Plant II Binding Plate. Do not moisten the rims of the individual wells while dispensing the samples.

Optional: Seal openings of the binding plate with a Gas-permeable Foil.

6 Bind DNA to silica membrane

Place the NucleoSpin® Plant II Binding Plate stacked on an MN Square-well Block in the rotor buckets. Centrifuge at **5,600–6,000 x g** for 5 min.

Typically, lysates will pass through the columns within 1 min. The centrifugation process can be extended to 20 min, if the lysates have not passed completely.

7 Wash silica membrane

1st wash

Add **400 µL PW1** to each well of the NucleoSpin® Plant II Binding Plate.

Seal plate with a Gas-permeable Foil and centrifuge again at **5,600–6,000 x g** for 2 min. Place NucleoSpin® Plant II Binding Plate on a new MN Square-well Block.

2nd wash

Add **700 µL PW2** to each well of the NucleoSpin® Plant II Binding Plate.

Seal plate with a Gas-permeable Foil and centrifuge again at **5,600–6,000 x g** for 2 min.

3rd wash

Add **700 µL PW2** to each well of the NucleoSpin® Plant II Binding Plate.

Seal plate with a Gas-permeable Foil and centrifuge again at **5,600–6,000 x g** for 10 min. Place NucleoSpin® Plant II Binding Plate on a new MN Square-well Block.

Note: For critical ethanol-sensitive applications, it is recommended to prolong the centrifugation time up to 15 min or incubate at higher temperature. Remove the adhesive foil and place the NucleoSpin® Plant II Binding Plate into an incubator for 20 min at 37 °C to evaporate residual ethanol.

8 Elute DNA

Place NucleoSpin® Plant II Binding Plate on the Rack of Tube Strips. Dispense **100 µL preheated Buffer PE (70 °C)** to each well of the NucleoSpin® Plant II Binding Plate. Dispense the buffer directly onto the membrane.

Optional: Incubate for 2 min at 70 °C before centrifugation.

Centrifuge at **5,600–6,000 x g** for 2 min. Remove the NucleoSpin® Plant II Binding Plate from the Rack of Tube Strips.

For optimal yield it is recommended to repeat this step once (incubation of Buffer PE on the membrane not required).

Yields will be 10–20 % higher when eluting with 2 x 100 µL Buffer PE depending on the total amount of DNA. However, the concentration of DNA will be much lower than with 100 µL.

Note: Elution can be done with TE buffer (at least pH 8.0) as well. Elution efficiency will decrease when using elution buffers with pH ≤ 8.0.

5.2 NucleoSpin® 96 Plant II – vacuum processing

- For hardware requirements, refer to section 2.3.
- For detailed information regarding the vacuum manifold setup, see page 23.
- For detailed information on each step, see page 23.
- For use of the NucleoSpin® 96 Plant II Core Kit (REF 740468.4), refer to section 2.4 regarding recommended accessories.

Before starting the preparation:

- Check if Buffer PW2 and RNase A were prepared according to section 3.
- Set incubator or oven to 65 °C.
- Equilibrate Buffer PE to 70 °C.

Protocol at a glance


1	Homogenize samples	Up to 100 mg wet or 20 mg lyophilized plant tissue 5,600–6,000 x g, 2 min
2a	Cell lysis using Buffer PL1	500 µL PL1 10 µL RNase A Mix 65 °C, 30 min Proceed with step 3
2b	Cell lysis using Buffer PL2 and PL3	400 µL PL2 Mix 65 °C, 30 min 100 µL PL3 Mix and incubate on ice for 5 min Proceed with step 3
3	Clear lysate by centrifugation	5,000–6,000 x g 20 min
4	Adjust binding conditions	Mix 450 µL PC with 400 µL cleared lysate
5	Transfer lysate to NucleoSpin® Plant II Binding Plate	

6	Bind DNA to silica membrane of the NucleoSpin® Plant II Binding Plate	-0.2 to -0.4 bar* (2 min)
7	Wash and dry silica membrane	400 µL PW1 700 µL PW2 700 µL PW2–0.4 bar* (1 min each step)
		Remove MN Wash Plate
		Dry silica membrane (10 min, maximum vacuum)
8	Elute DNA	100 µL PE (incubate 2 min) -0.4 bar* (2 min) Repeat once


* Reduction of atmospheric pressure

Setup of vacuum manifold:

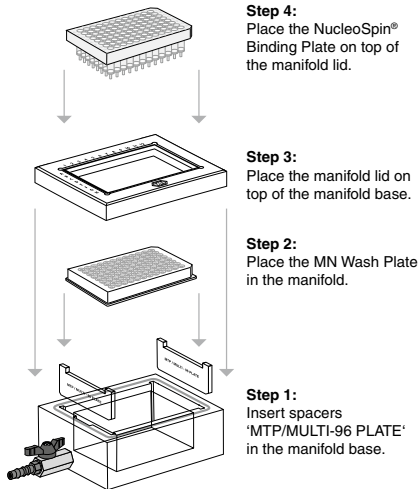Binding / Washing steps



Elution step

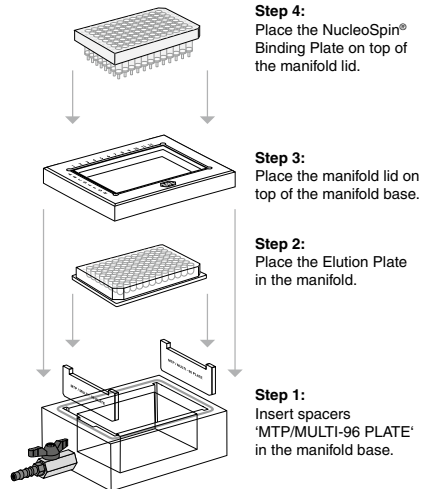

Detailed protocol

- For hardware requirements, refer to section 2.3.
- For detailed information regarding the vacuum manifold setup, see page 23.
- For use of the NucleoSpin® 96 Plant II Core Kit (REF 740468.4), refer to section 2.4 regarding recommended accessories.

Before starting the preparation:

- Check if Buffer PW2 and RNase A were prepared according to section 3.
- Set incubator or oven to 65 °C.
- Equilibrate Buffer PE to 70 °C.
- For detailed information regarding the vacuum manifold set-up see page 23.

1 Homogenize samples

Fill **up to 100 mg wet plant tissue** (or up to **20 mg dried, for example lyophilized, plant tissue**) into each tube of the Tube Strips. Add one 3 mm diameter steel bead to each tube. Close the tubes with Cap Strips. Freeze samples in liquid nitrogen. Disrupt cells by vigorous shaking using a mixer mill. Spin at **5,600 x g** for **2 min** and remove Cap Strips.

For further processing use either Buffer PL1 (2 a) or Buffers PL2/PL3 (2 b)!

2a Cell lysis using Buffer PL1

Add **500 µL Buffer PL1** and **10 µL RNase A** to each sample. Close tubes again using new Cap Strips (supplied). Mix by vigorous **shaking** for **15–30 s**. **Spin briefly** for **30 s** at **1,500 x g** to collect any sample from the Cap Strips. Incubate samples at **65 °C** for **30 min**.

Depending on plant sample and available methods, Buffer PL1 and RNase A may be added to the plant material before homogenization by the appropriate mechanical method.

Proceed with step 3.

2b Cell lysis using Buffer PL2 and PL3

Add **400 µL Buffer PL2** and **10 µL RNase A** to each sample. Close tubes again using new Cap Strips (supplied). Mix by vigorous **shaking** for **15–30 s**. **Spin briefly** for **30 s** at **1,500 x g** to collect any sample from the Cap Strips. Incubate samples at **65 °C** for **30 min**.

Depending on plant sample and available methods, Buffer PL2 and RNase A may be added to the plant material before homogenization by the appropriate mechanical method.

Open tubes, add **100 µL Buffer PL3**, close tubes, mix thoroughly, and incubate for 5 min on ice to precipitate SDS completely.

3 Clear lysate by centrifugation

Centrifuge the samples for **20 min** at **full speed (5,600–6,000 x g)**. Remove Cap Strips.

4 Adjust binding conditions

Predispense **450 µL Binding Buffer PC** to each well of a MN Square-well Block. Add **400 µL cleared lysate** of each sample and mix by repeated pipetting up and down. Mix at least 3 times.

5 Transfer lysate to NucleoSpin® Plant II Binding Plate

Place waste tray into manifold base. Insert spacers labeled 'MTP/MULTI-96 PLATE' notched side up into NucleoVac and place the MN Wash Plate on them. Close manifold and place NucleoSpin® Plant II Binding Plate on top of the manifold.

Transfer samples from the previous step into the wells of the NucleoSpin® Plant II Binding Plate. Do not moisten the rims of the individual wells while dispensing the samples.

6 Bind DNA to silica membrane

Apply vacuum of **-0.2 to -0.4 bar*** to allow samples to pass through the membrane. Flowthrough rate should be about 1–2 drops per second. Adjust vacuum strength accordingly. Finally, release the vacuum.

7 Wash silica membrane

1st wash

Add **400 µL PW1** to each well of the NucleoSpin® Plant II Binding Plate and apply vacuum of **-0.2 to -0.4 bar*** until the buffer has passed the membrane completely. Release the vacuum.

2nd wash

Add **700 µL PW2** to each well of the NucleoSpin® Plant II Binding Plate and apply vacuum of **-0.2 to -0.4 bar*** until the buffer has passed the membrane completely. Release the vacuum.

3rd wash

Add **700 µL PW2** to each well of the NucleoSpin® Plant II Binding Plate and apply vacuum of **-0.2 to -0.4 bar*** until the buffer has passed the membrane completely. Release the vacuum.

Remove MN Wash Plate and waste tray.

Reassemble the vacuum manifold and dry the membrane by applying maximum vacuum (**-0.6 bar***) for 10 minutes.

* Reduction of atmospheric pressure

8 Elute DNA

Insert spacers 'MICROTUBE RACK' into the vacuum manifold base. Place the Rack of Tube Strips into the manifold base. Close the manifold and insert the NucleoSpin® Plant II Binding Plate into the manifold top. Dispense **100 µL preheated Buffer PE (70 °C)** to each well of the NucleoSpin® Plant II Binding Plate. Dispense the buffer directly onto the membrane. Incubate at room temperature for 2 min. Apply vacuum of **-0.4 bar*** until the elution buffer has passed the membrane completely.

For optimal yield it is recommended to repeat this step once (incubation of Buffer PE on the membrane not required)

Yields will be 10–20 % higher when eluting with 2 x 100 µL Buffer PE depending on the total amount of DNA. However, the concentration of DNA will be much lower than with 100 µL.

***Note:** Elution can be done with TE buffer (at least pH 8.0) as well. Elution efficiency will decrease when using elution buffers with pH ≤ 8.0.*

* Reduction of atmospheric pressure

6 Appendix

6.1 Troubleshooting

Problem	Possible cause and suggestions
	<i>Homogenization of plant material was not sufficient</i> - For most species we recommend grinding with steel beads. Homogenization should be done thoroughly until the plant material is ground to a fine powder. In most cases this can be achieved by vigorous shaking for 3 x 60 s with occasional freezing in liquid nitrogen. - This problem can also be avoided by lyophilizing the material. This way, it will be easier to grind the material. <i>Extraction of DNA from plant material during lysis was not sufficient</i> - To obtain higher yields of DNA, the incubation time in lysis buffer can be prolonged (up to overnight). <i>Suboptimal lysis buffer was used</i> DNA yield is low - Lysis efficiencies of Buffer PL1 (CTAB) and Buffer PL2 (SDS) are different and depend on the plant species. Try both buffers in a side-by-side purification to find the best detergent system to lyse your plant material. <i>Sample contains too much RNA</i> - Add 10 µL of RNase A solution to the Lysis Buffer PL1 or PL2 before heat incubation. If this is not successful, add the enzyme to the cleared supernatant of step 3 and incubate for 30 min at 60 °C. <i>Sub-optimal Elution</i> - The DNA can be either eluted in higher volumes (up to 300 µL) or by repeating the elution step up to three times. Remember that the elution buffer must be preheated to 70 °C prior to elution. - Also check the pH of the elution buffer used, which should be in a range of pH 8–8.5. To ensure correct pH, use supplied elution Buffer PE.
DNA is degraded	<i>Sample was contaminated with DNase</i> Check bench, pipettes and storage of sample in order to avoid DNase contamination.

Problem	Possible cause and suggestions
DNA purity is low	<i>Sample contains DNA-degrading contaminants (e.g., phenolic compounds, secondary metabolites)</i>
	• Repeat washing step with Buffer PW1.
	<i>Elution buffer contains EDTA</i> • EDTA can disturb subsequent reactions. Use of water or supplied Elution Buffer PE is highly recommended.

6.2 Ordering information

Product	REF	Pack of
NucleoSpin® 96 Plant II	740663.2	1 x 96 preps
	740663.4	4 x 96 preps
	740663.24	24 x 96 preps
NucleoSpin® 96 Plant II Core Kit	740468.4	4 x 96 preps
NucleoSpin® 8 Plant II	740669	12 x 8 preps
	740669.5	60 x 8 preps
NucleoSpin® 8 Plant II Core Kit	740467.4	48 x 8 preps
MN Bead Plate Type D	740853.1	pack of 1
	740853.4	pack of 4
	740853.24	pack of 24
MN Bead Tube Type G	740817.50	50 pieces
Buffer PL1	740918	125 mL
Buffer Set PL2 / PL3 (100 mL Buffer PL2 + 25 mL Buffer PL3)	740919	1 set
Buffer PC	740937	125 mL
Buffer PW1	740938	125 mL
Buffer PW2 Concentrate (for 250 mL Buffer PW2)	740939	50 mL
RNase A (lyophilized)	740505	100 mg
	740505.50	50 mg
Proteinase K	740506	100 mg
MN Square-well Block	740476	4
	740476.24	24
MN Wash Plate	740479	4
	740479.24	24
Rack of Tube Strips (1 set consists of 1 rack, 12 strips with 8 tubes each, and 12 Cap Strips)	740477	4 sets
	740477.24	24 sets
Cap Strips	740478	48
	740478.24	288
Gas-permeable Foil	740675	50
Self adhering Foil	740676	50

Product	REF	Pack of
NucleoVac 96 Vacuum Manifold	740681	1
NucleoVac Vacuum Regulator	740641	1

Visit www.mn-net.com for more detailed product information.

6.3 Product use restriction / warranty

NucleoSpin® 96 Plant II (Core Kit) components are intended, developed, designed, and sold FOR RESEARCH PURPOSES ONLY, except, however, any other function of the product being expressly described in original MACHEREY-NAGEL product leaflets.

MACHEREY-NAGEL products are intended for GENERAL LABORATORY USE ONLY! MACHEREY-NAGEL products are suited for QUALIFIED PERSONNEL ONLY! MACHEREY-NAGEL products shall in any event only be used wearing adequate PROTECTIVE CLOTHING. For detailed information please refer to the respective Material Safety Data Sheet of the product! MACHEREY-NAGEL products shall exclusively be used in an ADEQUATE TEST ENVIRONMENT. MACHEREY-NAGEL does not assume any responsibility for damages due to improper application of our products in other fields of application. Application on the human body is STRICTLY FORBIDDEN. The respective user is liable for any and all damages resulting from such application.

DNA/RNA/PROTEIN purification products of MACHEREY-NAGEL are suitable for *IN VITRO*-USES ONLY!

ONLY MACHEREY-NAGEL products specially labeled as IVD are also suitable for *IN VITRO*-diagnostic use. Please pay attention to the package of the product. *IN VITRO*-diagnostic products are expressly marked as IVD on the packaging.

IF THERE IS NO IVD SIGN, THE PRODUCT SHALL NOT BE SUITABLE FOR *IN VITRO*-DIAGNOSTIC USE!

ALL OTHER PRODUCTS NOT LABELED AS IVD ARE NOT SUITED FOR ANY CLINICAL USE (INCLUDING, BUT NOT LIMITED TO DIAGNOSTIC, THERAPEUTIC AND/OR PROGNOSTIC USE).

No claim or representations is intended for its use to identify any specific organism or for clinical use (included, but not limited to diagnostic, prognostic, therapeutic, or blood banking). It is rather in the responsibility of the user or - in any case of resale of the products - in the responsibility of the reseller to inspect and assure the use of the DNA/RNA/protein purification products of MACHEREY-NAGEL for a well-defined and specific application.

MACHEREY-NAGEL shall only be responsible for the product specifications and the performance range of MN products according to the specifications of in-house quality control, product documentation and marketing material.

This MACHEREY-NAGEL product is shipped with documentation stating specifications and other technical information. MACHEREY-NAGEL warrants to meet the stated specifications. MACHEREY-NAGEL's sole obligation and the customer's sole remedy is limited to replacement of products free of charge in the event products fail to perform as warranted. Supplementary reference is made to the general business terms and conditions of MACHEREY-NAGEL, which are printed on the price list. Please contact us if you wish to get an extra copy.

There is no warranty for and MACHEREY-NAGEL is not liable for damages or defects arising in shipping and handling (transport insurance for customers excluded), or out of accident or improper or abnormal use of this product; defects in products or components not manufactured by MACHEREY-NAGEL, or damages resulting from such non-MACHEREY-NAGEL components or products.

MACHEREY-NAGEL makes no other warranty of any kind whatsoever, and SPECIFICALLY DISCLAIMS AND EXCLUDES ALL OTHER WARRANTIES OF ANY KIND OR NATURE WHATSOEVER, DIRECTLY OR INDIRECTLY, EXPRESS OR IMPLIED, INCLUDING, WITHOUT LIMITATION, AS TO THE SUITABILITY, REPRODUCTIVITY, DURABILITY, FITNESS FOR A PARTICULAR PURPOSE OR USE, MERCHANTABILITY, CONDITION, OR ANY OTHER MATTER WITH RESPECT TO MACHEREY-NAGEL PRODUCTS.

In no event shall MACHEREY-NAGEL be liable for claims for any other damages, whether direct, indirect, incidental, compensatory, foreseeable, consequential, or special (including but not limited to loss of use, revenue or profit), whether based upon warranty, contract, tort (including negligence) or strict liability arising in connection with the sale or the failure of MACHEREY-NAGEL products to perform in accordance with the stated specifications. This warranty is exclusive and MACHEREY-NAGEL makes no other warranty expressed or implied.

The warranty provided herein and the data, specifications and descriptions of this MACHEREY-NAGEL product appearing in MACHEREY-NAGEL published catalogues and product literature are MACHEREY-NAGEL's sole representations concerning the product and warranty. No other statements or representations, written or oral, by MACHEREY-NAGEL's employees, agent or representatives, except written statements signed by a duly authorized officer of MACHEREY-NAGEL are authorized; they should not be relied upon by the customer and are not a part of the contract of sale or of this warranty.

Product claims are subject to change. Therefore please contact our Technical Service Team for the most up-to-date information on MACHEREY-NAGEL products. You may also contact your local distributor for general scientific information. Applications mentioned in MACHEREY-NAGEL literature are provided for informational purposes only. MACHEREY-NAGEL does not warrant that all applications have been tested in MACHEREY-NAGEL laboratories using MACHEREY-NAGEL products. MACHEREY-NAGEL does not warrant the correctness of any of those applications.

Last updated: 07 / 2010, Rev. 03

Please contact:

MACHEREY-NAGEL GmbH & Co. KG

Tel.: +49 (0) 24 21 969 270

e-mail: tech-bio@mn-net.com

Trademarks / disclaimer:

NucleoSpin® is a registered trademark of MACHEREY-NAGEL GmbH & Co. KG

All used names and denotations can be brands, trademarks, or registered labels of their respective owner – also if they are not special denotation. To mention products and brands is only a kind of information (i.e., it does not offend against trademarks and brands and can not be seen as a kind of recommendation or assessment). Regarding these products or services we can not grant any guarantees regarding selection, efficiency, or operation.



www.mn-net.com

MACHEREY-NAGEL



MACHEREY-NAGEL GmbH & Co. KG · Neumann-Neander-Str. 6-8 · 52355 Düren · Germany

DE/International:

Tel.: +49 24 21 969-0
Fax: +49 24 21 969-199
E-mail: info@mn-net.com

CH:

Tel.: +41 62 388 55 00
Fax: +41 62 388 55 05
E-mail: sales-ch@mn-net.com

FR:

Tel.: +33 388 68 22 68
Fax: +33 388 51 76 88
E-mail: sales-fr@mn-net.com

US:

Tel.: +1 484 821 0984
Fax: +1 484 821 1272
E-mail: sales-us@mn-net.com



A035556/0800.5